

Specimen Holders



**Specimen holders for
unmounted or
mounted specimens
of all shapes and sizes**



- Specimen holders for all sample shapes and sizes
- Specimen holders for mounted, unmounted or irregular specimens
- Interchangeable from one machine to another
- Holders available in 140, 160 and 200 mm dia.
- Stainless steel or light-weight aluminium holders
- Dual clamping for increased capacity
- Up to 12 specimens can be clamped at a time
- Adjustable hole diameter provides increased flexibility
- Save valuable time by easy and fast clamping





The basis of a perfect preparation

Automatic preparation of metallographic specimens requires a level clamping of the specimens in specimen holders.

Struers offers a wide range of specimen holders for clamping of 3 or more mounted or unmounted specimens. All kinds of shapes and sizes can be clamped.

All specimen holders have the same type of coupling, which makes them interchangeable from one Struers machine to another.

The range of specimen holders includes 140, 160 and 200 mm dia. sizes in both stainless steel and aluminium.



The specimen holder with rectangular holes can hold specimens up to 70x53 mm



All kinds of shapes and sizes can be clamped in Struers specimen holders

Holders for mounted or unmounted specimens


For mounted, unmounted or irregular specimens, Struers offers a number of special specimen holders with teardrop or rectangular holes.

Specimen holders with rectangular holes can hold specimens up to 70x53 mm and you avoid having to cut the specimen into smaller pieces.



Specimen holder with clamps

This holder comes in handy when you need to clamp wire ends or plates. Three individual, detachable clamps are fitted on a specimen holder.

This allows for easier fixing of several specimens, preparation of both faces, and even cutting of specimens fitted in the clamps prior to preparation.

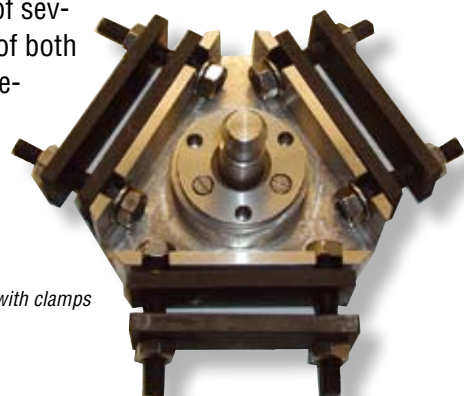

Specimen holder with clamps

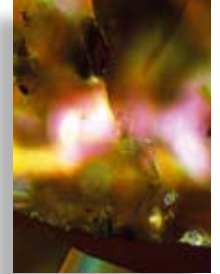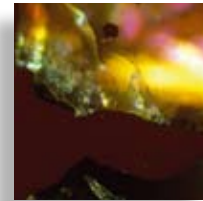

Tailor-made solutions

Plane, blank specimen holders for adhesion of typically large plane parallel specimens are also available.

These specimen holders are also used to produce tailor-made solutions for customers with special requirements.

Application

For high capacity machines like MAPS or AbraPol machines, holders with a diameter up to 200 mm can be used.

On the TegraSystem, holders with a diameter up to 160 mm (200 mm) can be used on TegraForce-3 and TegraForce-5.

Specimen mover plates

Struers offers a wide range of specimen mover plates for preparation of single specimens.

Dual-Clamping - saves time and increases your capacity

One screw fastens 2 specimens

Struers programme also includes light-weight aluminium specimen holders designed with a dual-clamping principle: One screw fastens two specimens. This feature facilitates faster clamping and increases the capacity of the specimen holder.

Dual-clamping also allows a certain tolerance for the specimen diameter, which makes these specimen holders well suited for all specimens.

The specimen holders are available in three different sizes, 140, 160 and 200 mm dia., and can hold up to 12 specimens at a time.



Dual-clamping - two specimens are fastened with a single screw

Inserts - increased flexibility

Inserts of reinforced nylon allow you to adjust the hole diameter easily. You just click the insert into the holder and it is ready to be used for smaller specimens. This is especially useful if you only prepare smaller specimens occasionally.

This unique feature allows you to place specimens of different sizes in the same specimen holder, thus adding to the flexibility and reducing the amount of specimen holders needed to do the job.

Three different inserts are available to change the hole size from:

50 mm (2") to 40 mm (1½") dia.



40 mm (1½") to 30 mm (1¼") dia.

30 mm (1¼") to 25 mm (1") dia.

Uniforce - for fast and accurate levelling

Efficient clamping of the specimens saves valuable time in the following preparation steps.

With Uniforce Struers offers fast and accurate levelling of specimens in 140, 160 and 200 mm dia. specimen holders. Precise levelling will reduce your plane grinding time to a minimum and ensure the best possible conditions for quick metallographic preparation of high quality.

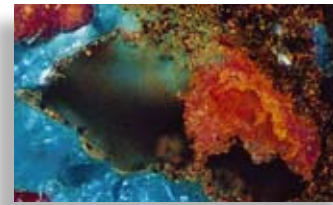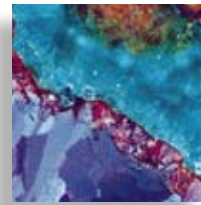

Inserts provide easy adjustment of the hole diameter



With Uniforce, the projection of specimens from the holder can be adjusted to 3 or 6 mm, with 3 mm being the standard.

For more detailed information on Uniforce, ask for a separate brochure.



Uniforce

Specimen Holders

Drawing	Dia.	Specimens	Size	Code
	140 mm	8	24-26 mm (1") dia.	SAXCA
		6	29-32 mm (1 1/4") dia.	SAXCE
			30 mm (1 1/4") to 25 mm (1") dia.	SAXIT
		6	10-32 mm dia.	PEDST
				PEDAL
		4	10-40 mm dia.	PEDIS
				PEDFI
				MAXAD
	160 mm	10	24-26 mm (1") dia.	SAXDA
		8	29-32 mm (1 1/4") dia.	SAXDE
		6	38-40 mm (1 1/2") dia.	SAXDI
			40 mm (1 1/2") to 30 mm (1 1/4") dia.	SAXIO
			30 mm (1 1/4") to 25 mm (1") dia.	SAXIT
		6	10-40 mm dia.	MAXCY
				MAXLU
		12	10-25.5 mm dia.	MAXDU
		6	25x34 mm	MAXYK
		3	48x40 mm *	MAXDI
		3	24x52 mm	MAXAB
				MAXSO
				MAXAS

Drawing	Dia.	Specimens	Size	Code
	200 mm	12	29-32 mm (1 1/4") dia.	SAXFE
		8	38-40 mm (1 1/2") dia.	SAXFI
		6	49-52 mm (2") dia.	SAXFO
			50 mm (2") to 40 mm (1 1/2") dia.	SAXIF
			40 mm (1 1/2") to 30 mm (1 1/4") dia.	SAXIO
			30 mm (1 1/4") to 25 mm (1") dia.	SAXIT
		5	27-55 mm dia.	MAXUK
			3 specimens of max. 70 mm dia.	MAXRU
			1 specimen of max: outer dia.: 180 mm inner dia.: 40 mm	MAXOT
		3	53x70 mm **	MAXOT
				MAXON
				MAXAN

Code symbols	Legend
Code	Reinforced nylon
Code	Stainless steel
Code	Aluminium

* The MAXDI specimen holder can also be used for specimens mounted in UnoForm (UNOQA)

** The MAXOT specimen holder can also be used for specimens mounted in UnoForm (UNOQI)



Struers' products are subject to constant product development. Therefore, we reserve the right to introduce changes in our products without notice.



Struers A/S
Pederstrupvej 84
DK-2750 Ballerup, Denmark
Phone +45 44 600 800
Fax +45 44 600 801
struers@struers.dk

USA
Struers Inc.
24766 Detroit Road
Westlake, OH 44145-1598
Phone +1 440 871 0071
Fax +1 440 871 8188
info@struers.com

CANADA
Struers Ltd.
7275 West Credit Avenue
Mississauga, Ontario L5N 5M9
Phone +1 905-814-8855
Fax +1 905-814-1440
info@struers.com

SWEDEN
Struers A/S
Ekbacksvägen 22, 3 tr
P.O. Box 11085
SE-168 69 Bromma
Telefon +46 (0)8 447 53 90
Telefax +46 (0)8 447 53 99
struers@struers.dk

FRANCE
Struers S.A.S.
370, rue du Marché Rollay
F- 94507 Champigny
sur Marne Cedex
Téléphone +33 1 5509 1430
Télécopie +33 1 5509 1449
struers@struers.fr

NETHERLAND/BELGIE
Struers GmbH Nederland
Electraweg 5
NL-3144 CB Maassluis
Tel. +31 (0) 10 599 72 09
Fax +31 (0) 10 599 72 01
glen.van.vugt@struers.de

BELGIE (Wallonie)
Struers S.A.S.
370, rue du Marché Rollay
F- 94507 Champigny
sur Marne Cedex
Téléphone +33 1 5509 1430
Télécopie +33 1 5509 1449
struers@struers.fr

UNITED KINGDOM
Struers Ltd.
Unit 25a
Monkspath Business Park
Solihull
B90 4NZ
Phone +44 0121 745 8200
Fax +44 0121 733 6450
info@struers.co.uk

IRELAND
Struers Ltd.
Unit 25a
Monkspath Business Park
Solihull
B90 4NZ
Phone +44 (0)121 745 8200
Fax +44 (0)121 733 6450
info@struers.co.uk

JAPAN
Marumoto Struers K.K.
Takara 3rd Building
18-6, Higashi Ueno 1-chome
Taito-ku, Tokyo 110-0015
Phone +81 3 5688 2914
Fax +81 3 5688 2927
struers@struers.co.jp

CHINA
Struers Ltd.
Office 702 Hi-Shanghai
No. 970 Dalian Road
Shanghai 200092, P.R. China
Phone +86 (21) 5228 8811
Fax +86 (21) 5228 8821
struers.cn@struers.dk

DEUTSCHLAND
Struers GmbH
Karl-Arnold-Strasse 13 B
D-47877 Willich
Telefon +49(0)2154 486-0
Telefax +49(0)2154 486-222
verkauf.struers@struers.de

ÖSTERREICH
Struers GmbH
Zweigniederlassung Österreich
Ginzkeyplatz 10
A-5020 Salzburg
Telefon +43 662 625 711
Telefax +43 662 625 711 78
stefan.lintschinger@struers.de

SCHWEIZ
Struers GmbH
Zweigniederlassung Schweiz
Weissenbrunnstrasse 41
CH-8903 Birmensdorf
Telefon +41 44 777 63 07
Telefax +41 44 777 63 09
rudolf.weber@struers.de

CZECH REPUBLIC
Struers GmbH
Organizační složka
Havlíčková 361
CZ-252 63 Roztoky u Prahy
Tel: +420 233 312 625
Fax: +420 233 312 640
david.cernicky@struers.de

POLAND
Struers Sp. z o.o.
Oddział w Polsce
ul. Lirowa 27
PL-02-387 Warszawa
Tel. +48 22 824 52 80
Fax +48 22 882 06 43
grzegorz.uszynski@struers.de

HUNGARY
Struers GmbH
Magyarországi fióktelep
Puskás Tivadar u. 4
H-2040 Budaörs
Phone +36 (23) 428-742
Fax +36 (23) 428-741
zoltan.kiss@struers.de

SINGAPORE
Struers A/S
627A Aljunied Road,
#07-08 BizTech Centre
Singapore 389842
Phone +65 6299 2268
Fax +65 6299 2661
struers.sg@struers.dk

www.struers.com